



EPIDUROSCOPY


A new way to look at back pain



Via Martiri di Cefalonia, 28 - 63822 Porto San Giorgio (FM) - Italy
Tel. +39 0734 674379 Fax +39 0734 678487 www.alfamedsrl.com - info@alfamedsrl.com



Handle:
- Ergonomic
- Bi-directional guide mechanism which allows to move and direct the tip into the epidural space.

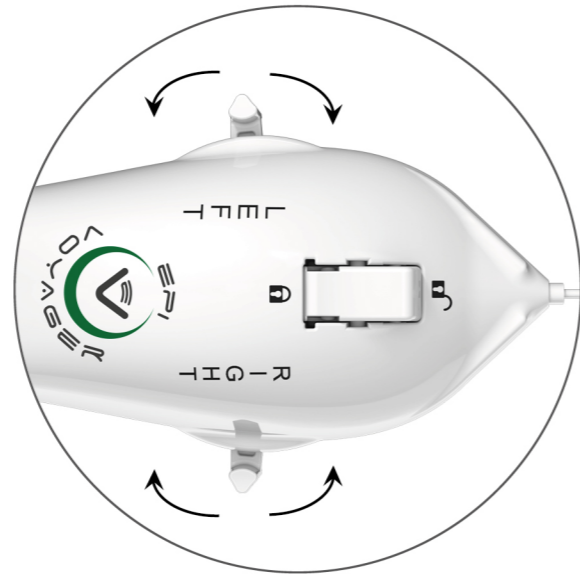

Infusion channels for simultaneous administration

4 WORKING CHANNELS

STYLET
Block valve to prevent displacement




EPI-VOYAGER

The steerable Videoguided EPI-VOYAGER Catheter is a device for spinal endoscopy to be used in the anterior and posterior epidural space of lumbo-sacral area. EPI-VOYAGER features 4 working channels.

Epidural Balloon Catheter



Tip infusion

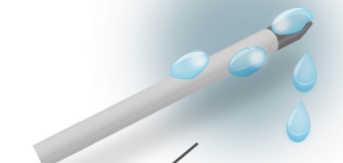

Flexible, atraumatic and radiopaque tip



Grasper

Flexible and Radiopaque catheter body with markers



Flexible steerable tip

REF	CATHETER LENGTH	CATHETER OUTER DIAMETER Ø	TOOL CHANNEL INNER DIAMETER Ø	OPTIC CHANNEL INNER DIAMETER Ø
EA3010	30cm	3.3 mm	10Fr	1.3 mm

"GAP" Epidural Balloon Catheter

Epidural Balloon Catheter, radiopaque, dedicated to adhesions removal. Certified for Epidural Use.

REF	DESCRIPTION
GA6035	GAP - Epidural Balloon Catheter - 60cm x 3,5Fr

